




3 Well Midi Crystallization Plate in UVXPO (Swissci) - 10/cs


Sitting drop crystallization

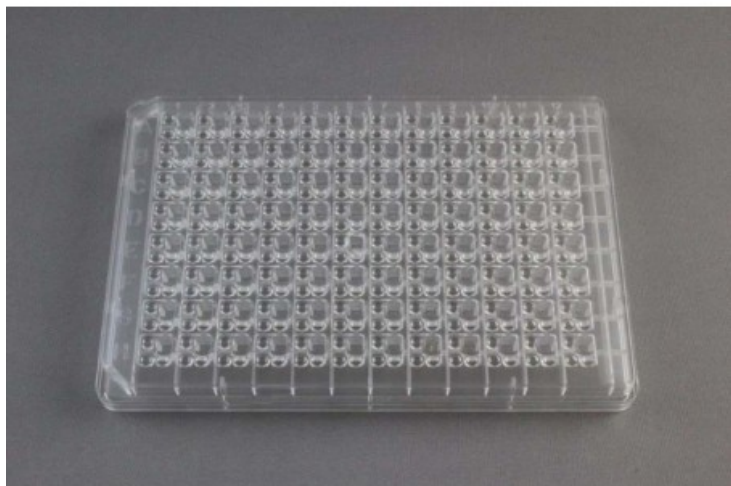

Product codes:

Reference: HR 3193

Product gallery:

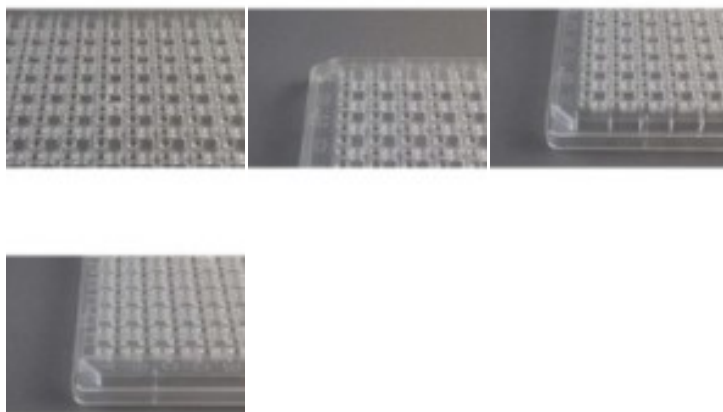

Product description:

- Lens effect drop wells for enhanced optics
- Drop volume: 50 nl to 5 μ l
- Micro-numbering alongside drop volumes
- Reservoir volume: 20 to 40 μ l
- ANSI/SLAS 1-2004 Standard
- Drop well allows easier crystal harvesting
- Three drop wells per reservoir
- Midi profile for easier storage
- UV compatible

The SWISSCI 3 Lens Crystallization Plate is an ideal product for protein crystallization presented in a 96-well plate format.

The SWISSCI UVXPO 3 Lens Crystallisation Plate is made from a new polymer mixture that means crystals can be viewed in UV, Polar and Visible light. UVXPO is short for UV compatible as well as excellent uniform (zero) background when using cross polarised light.



The SWISSCI 3 Lens Crystallisation Plate is an ideal product for protein crystallisation presented in a 96-well plate format.

The Advantages of the 3 Lens Crystallization Plate

Easy crystal retrieval

Raised wide wells make the crystal mounting especially easy.

Easy viewing

Crystals can be viewed in UV, Polar and Visible light.

Uniform (zero) background when using cross polarised light.

The wells are a wide conical shape and have a lens effect for perfect illumination.

The micro-numbering ensures you will never get lost again (visible by microscope).

The optically superior polymer is UV transmissible and may be used to differentiate between salt and protein crystals.

Better sealing

Wide partition walls between the wells give plenty of area for good sealing with tape.

No central bending occurs in this very robust structure.

Excellent long term storage - no sample evaporation.

Wide range of volumes

Typical volumes are 20-40 μ l of reservoir and 50 nl-5 μ l drop size. The 288 optical wells offer three times the number of experimental constructions.

ANSI/SLAS 1-2004 Standard

The plates are designed to the 96-well ANSI/SLAS 1-2004 standard for all common holders and external numbering (A-H, 1-12) with corner location make the plate easy to use in a robotic sampler.

The 3 Lens plate is suitable for centrifugation.

The unique 3 drop protein crystallisation plate offers a new way of sitting-drop crystallography.

The 288 wells are optically perfect designed to observe crystals under a microscope.

3 Lens plate features

Easy to fill 96-lens structure - typically 30 microlitres per buffer lens.

Standard ANSI/SLAS 1-2004 Standard form for automation. Optically perfect wells with lens effect for a better view with microscope. Microscopic identifier within wells, simplifying the orientation under the microscope. Maximum volume of the buffer reservoir is 40 microlitres, typically 10 nanolitres is used for the protein well filling. Grown crystals are easy to identify and to remove from lens due to a low-binding polymer. Plate with 3 lenses for each sample, better growing security with triplicates or the ability to use lens two and lens three as mixing stations. Lenses fill without micro-droplets jumping out due to static effects. Wells can be individually sealed with a perfectly flat upper surface - there is a large land area designed to ensure integrity of each lens section. The profile allows for easier storage. Low volume buffer well enables savings on reagents. Triplicates allow for 288 constructs per plate.



Shop name

Per maggiori informazioni visita il sito <https://hamptonresearch.com/>

Product features:

CRF - TIPO: 3 Well Midi UVXPO;